



Get Life-like Translucency With Lava™ from Thoele

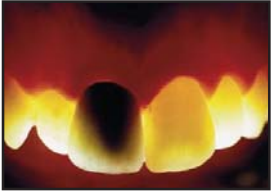

Of all the various all-ceramic crown and bridge materials out there, Lava™ from 3M continues to be the zirconia brand most asked for by dentists. Probably the most thoroughly researched and tested zirconia substrate on the market, Lava serves as the benchmark from which all other brands are compared. Lava is also the most translucent brand of zirconia, rendering life-like esthetics. Combine this with 1,100 MPa flexural strength and you can see why Thoele recommends Lava as an esthetic alternative to porcelain-fused-to-metal restorations. Please call us for more information.

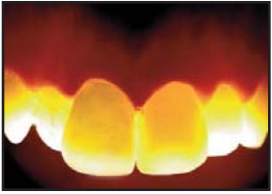

Call us about the Lava promotion from 3M and we'll send you coupons for \$75 off each of your next five Lava cases.

New! SR Phonares™ Denture Teeth from Ivoclar

Touted as denture teeth for implants, consider using Phonares™ teeth in your next denture case. Phonares is a new nano-hybrid composite (NHC)-based tooth which outshines existing denture teeth as far as esthetics and wear are concerned. The new tooth moulds are characterized by a unique surface texture which enhances the life-like appearance of the teeth. Another special feature of these moulds is the adjusted proximal surfaces which ensure an optimum level of esthetics. The translucency, opalescence and fluorescence make these teeth particularly attractive. The dentist may choose between youthful, universal and mature, as well as, between soft and bold tooth moulds. The SR Phonares anterior tooth line comprises 18 maxillary moulds and 4 mandibular moulds in an elaborate 4-layer design. They are available in the 16 Vita A-D shades and 4 Bleach shades. The posterior tooth lines are available in full contour to suit implant use as well as both standard 18-22 degree setups and lingualized occlusal schemes.

For more information on Phonares denture teeth call Jeff at 800.899.1115

Free SpeedCem™ with Your First e.max® Press Case

While quantities last, Thoele will include a free SpeedCem cement kit with first time e.max Press prescriptions. Ivoclar's new SpeedCem is a self-adhesive dual-curing universal dental resin cement. This universal resin cement has applications similar to 3M's Unicem. If you have any questions, would like a free sample of SpeedCem, or a cementation guide for the various restorations, please give Frank a call.

Reminder: Our custom shade-taking hours are between 7:30 & 11:30 a.m. and between 12:30 & 3:30 p.m. or by appointment

WANTED: Your e-mail address.

If you'd like us to have the option of communicating with you electronically, please send an e-mail to Stacy at sspoden@thoeledental.com or myself at fbender@thoeledental.com. We will not share your information with anyone outside of our lab without your approval. E-mail is a great way to ask questions about your case. You can even utilize digital photos or video clips of your cases to enhance communication. Please help us to better serve you and your patients.

Prescribe Custom Implant Abutments from Thoele



How many different implant systems does dentistry really need? And how do you gather enough parts, tools, and experience to make providing implants predictable and profitable? At last count there were 10 "major" implant systems available in the US and at least that many minor/private label systems. Within each system there are often multiple brands with a variety of materials, platforms, and dimensions.

Help! Although we can't eliminate all the choices, we can simplify the restorative process so you and your team can concentrate on patient care, and not on parts and pieces. Custom implant abutments such as Atlantis™ help make it possible to achieve highly esthetic and cost effective results. Custom computer designed and milled abutments are available for most implant systems, have predictable costs and provide precise interface to the implant platform for both biological health and screw stability. We choose only proven custom abutment manufacturers that are FDA cleared and hold the scientific, engineering and research credentials to ensure the highest quality and precision of the abutment that preserves the integrity of the implant itself. We do not use the lowest priced custom abutments because they lack these credentials.

Other Atlantis benefits include:

- Esthetic correction of offset or angled fixture placement
- Improved crown retention with ideal coronal preps
- Custom anatomical shape to mimic tissue contours
- Reduced precious metal weight of the final restoration
- Zirconia abutments allow for esthetic margins at tissue level making complete removal of excess cement easier

Try Metal Free Crowns Specifically Designed for the Posterior



Thoele now offers BruxZir™ Total Zirconia crowns for posterior restorations. BruxZir is virtually unbreakable and withstands severe parafunctional activity. BruxZir is a full contoured solid zirconia crown or bridge that is milled, stained and glazed. Because BruxZir has no porcelain veneer, it is more resistant to chipping, cracking or breaking in the mouth than other ceramic and porcelain-fused-to-metal options.

Despite their extraordinary strength, BruxZir crowns are adjustable chairside with a fine grit diamond bur using a light touch with water and air spray. This will keep the restoration cool and avoid micro-fractures while making adjustments. BruxZir provides your patients with a full gold crown or bridge alternative that has the same preparation requirements, is flat priced, lower cost, more esthetic and metal free.

Call us at 800.899.1115 for more information.

Upcoming Events:

Star of the North Meeting: April 22-24 Booth #541

Please stop by our booth for an informative chat, treats and entry into our annual drawing!

NTI/Obstructive Sleep Apnea Appliances:

- September 29, 2010 at Westside Learning and Events Center, Sartell
- September 30, 2010 at The Lodge, Brainard Lakes

Keep an eye out for other CE seminars being developed on such topics as all-ceramic options, chairside shade taking, implants and Locators.

